



***TomoTherapy®***  
ACCURAY®


# **TOMOHD™ TREATMENT SYSTEM**

Product Specifications


## 1. Complete Treatment System

### Dimensions


• Height	8 ft., 3.25 in. (252.1 cm)
• Width	9 ft., 2.50 in. (280.6 cm)
• Length	15 ft., 2.25 in. (463.0 cm)
• Weight	10,000 lbs. (4,536 kg)

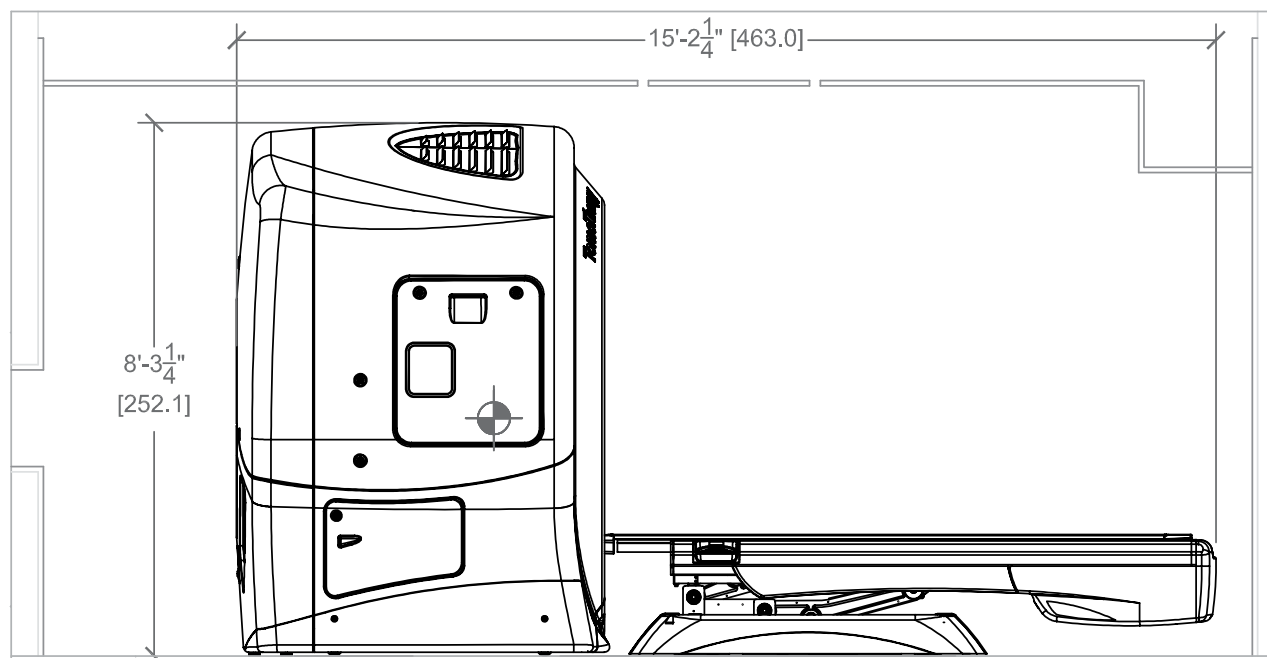

*Lateral View: System and minimum room dimensions*

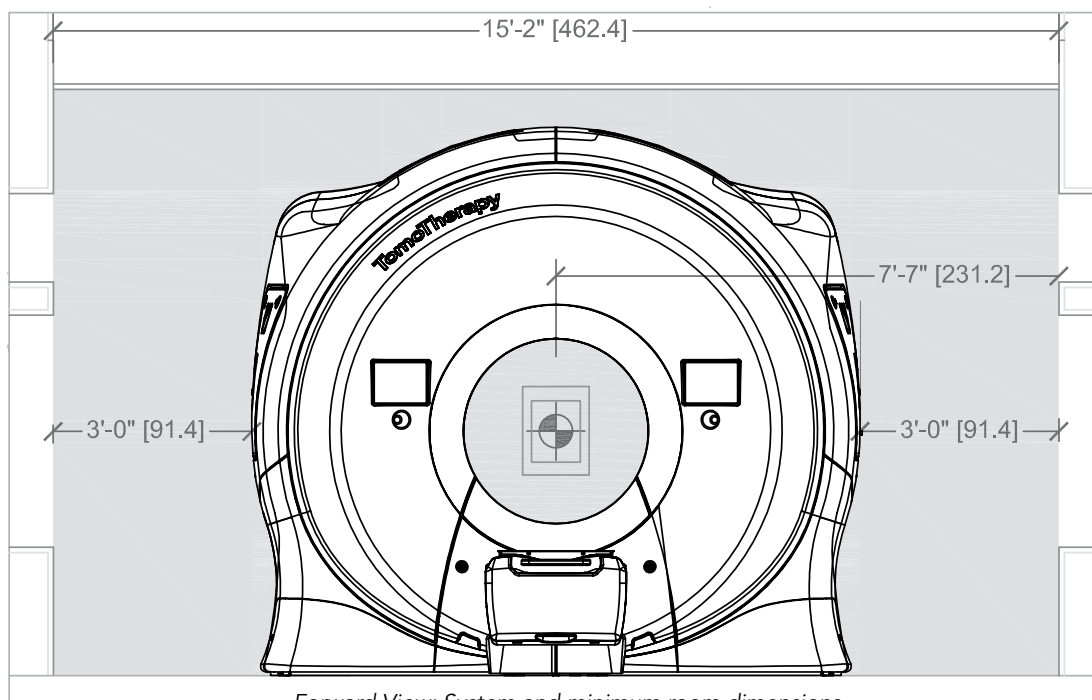

*Forward View: System and minimum room dimensions*

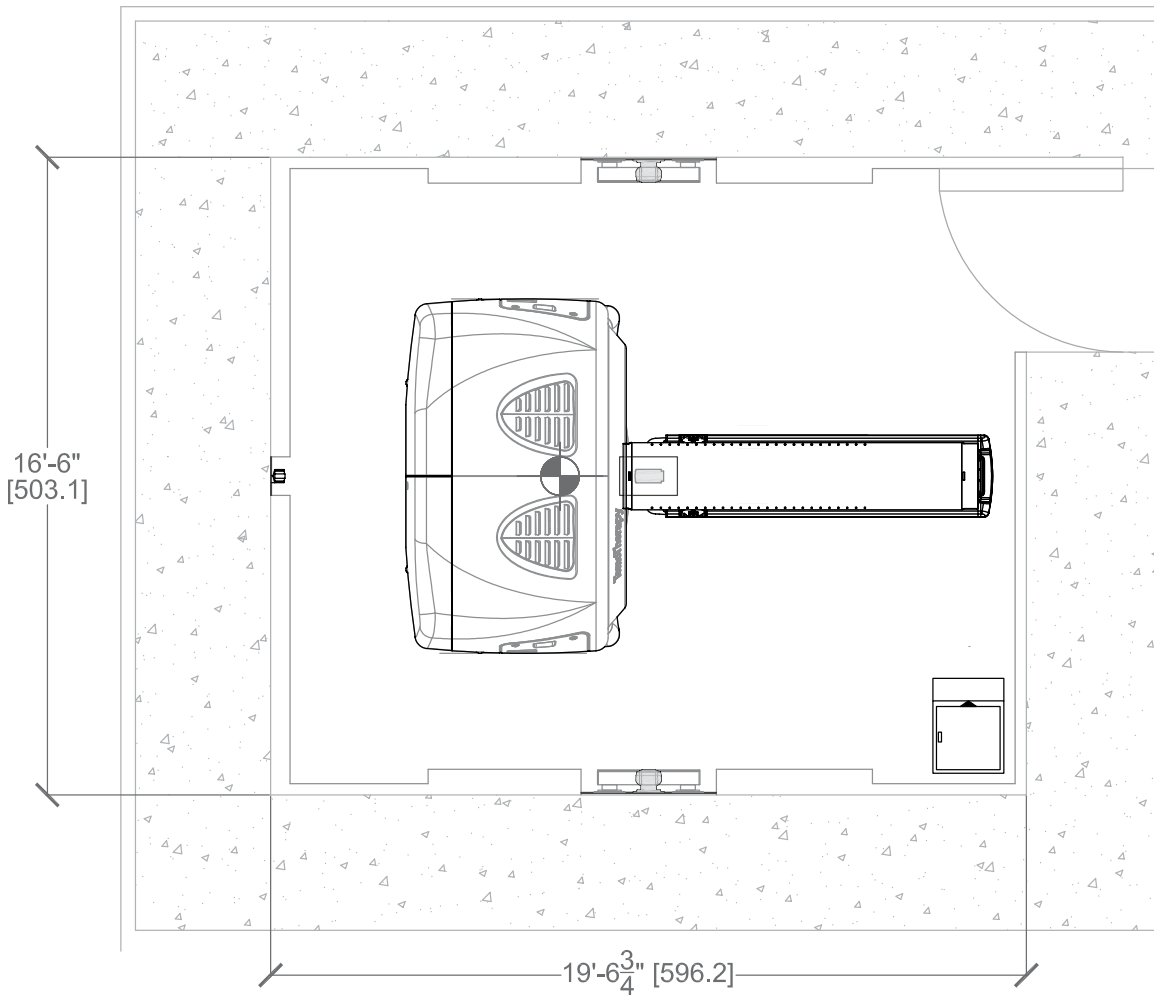

Top View: System and minimum room dimensions

## 2. Installation

### Minimum room size - requirements

- Height
- Width
- Length

8 ft., 10.13 in. (269.6 cm)

15 ft., 2.00 in. (462.4 cm)

19 ft., 6.75 in. (596.2 cm)

### Typical schedule (receipt of system to first treatment)

30 days, although completion has been achieved within two weeks

### Commissioning

Integrated system completely pre-commissioned in the factory and re-verified at customer site

## 3. Treatment Vault Environment

### SPECIFICATION

#### Line voltage, VAC

#### Ambient room temperature

#### Relative humidity

### REQUIREMENTS

480 volts, 3-phase; others may be supported

68 - 75 degrees Fahrenheit (20 - 24 degrees Celsius)

30% - 60%, noncondensing

## 4. Mechanical Features

### SPECIFICATION

#### Gantry

- Degrees of rotation
- Direction of rotation
- Rotational angle accuracy
- Speed of rotation
- Controls
- Source to axis distance
- Mechanical to radiation isocenter offset
- Mechanical isocenter stability
- Position indicators
- Isocenter height
- Cooling

### PERFORMANCE

Rotates around IEC-y axis, continuous rotation  
 Clockwise viewed from the foot of the couch  
 Within 0.5 degrees  
 Varies; dependent upon plan  
 Rotational speed set during treatment planning  
 85 cm  
 Included in beam model (within accuracy 0.25 mm)  
 < 0.4 mm  
 5-axis laser system  
 113 cm typical (dependent upon finished flooring)  
 Integrated onboard cooling system eliminates need for facility-chilled water loop

## 5. Photon Beam

### SPECIFICATION

#### Accelerator type

- Length

#### Microwave

- Power
- Source

#### Nominal dose rate at Dmax\*\*

#### Nominal Dmax\*\*

#### Percentage depth dose at 10 cm\*\*

#### Number of beamlets

#### Nominal energy

#### Nominal spot size

#### Penumbra of treatment beam - IEC-y axis

#### Field size range at isocenter

#### Maximum radiation field length

#### Treatment volume - TomoHelical™

#### Treatment volume - TomoDirect™

### PERFORMANCE

Standing Wave  
 0.3 meters

2.5 MW

Magnetron

850 cGy/min

1.5 cm

Typical 61.4%

Plan dependent; system capable of efficiently delivering tens of thousands of beamlets which allows for very high fidelity intensity modulation

6 MV, single energy

2 mm

The 80% - 20% penumbra widths are approximately 4.9 mm, 4.6 mm, and 4.1 mm (respectively for the 5.0 cm, 2.5 cm, and 1.0 cm field widths).

Selectable 1.0 cm x 0.625 cm to 1.0 cm x 40 cm  
 2.5 cm x 0.625 cm to 2.5 cm x 40 cm  
 5.0 cm x 0.625 cm to 5.0 cm x 40 cm


150 cm with Couch at height of isocenter plane

80 cm (transverse diameter) x 135 cm (longitudinal) for typical patient set-up; Actual treatment volume length is variable depending upon couch height; Note: Region of treatment used is determined by the planning CT image field of view (FOV) Provided the FOV is 80 cm or less and all patient anatomy is present in the planning image, the TomoHD system can import the image, plan, and treat.

40 cm (transverse diameter) x 135 cm (longitudinal) for typical patient set-up; Actual treatment volume length is variable depending upon couch height.

\*\* Performance based on 5 cm field size at SSD=85 cm.

Linear Accelerator (Linac)

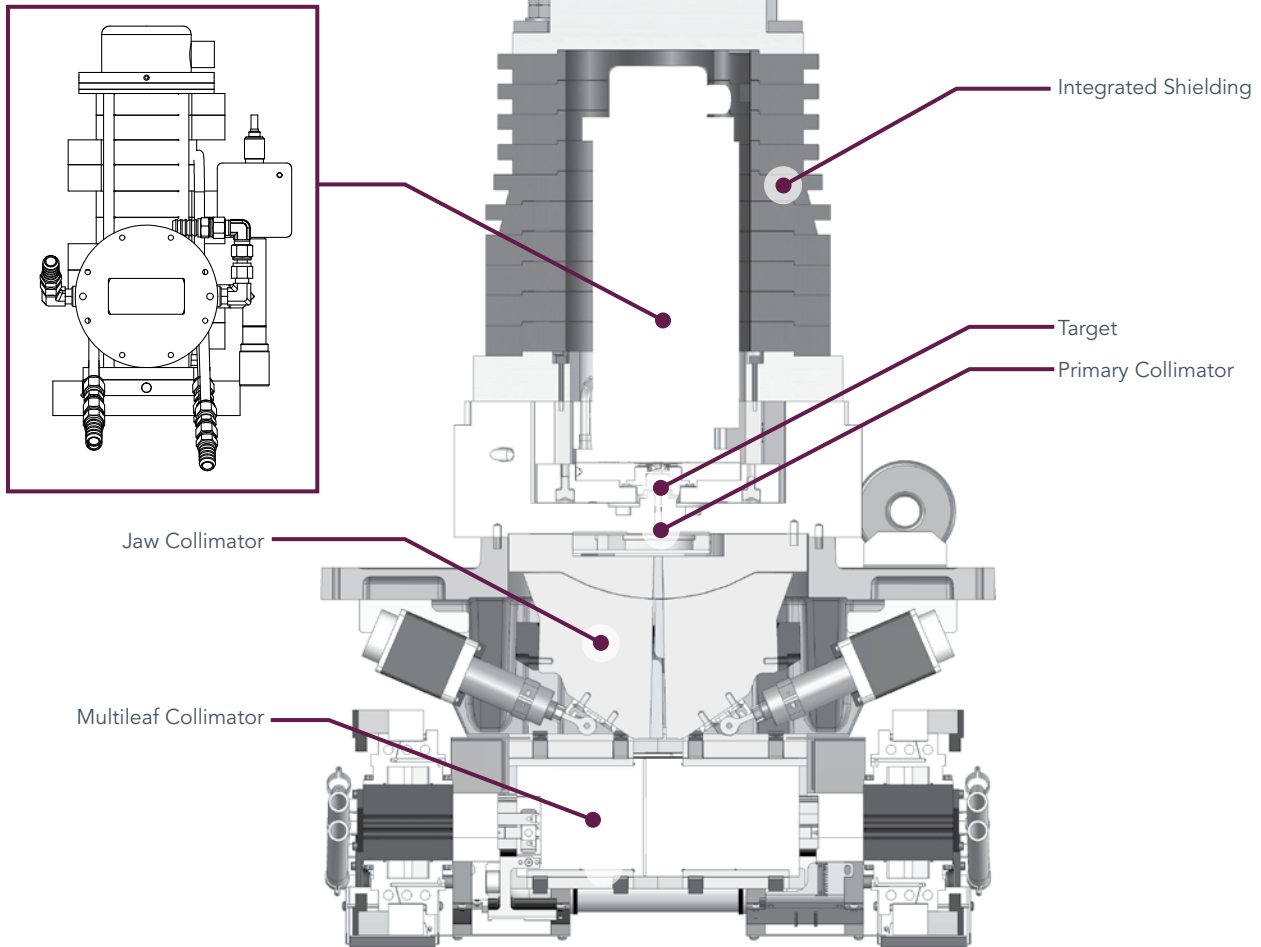

Accelerator and collimator subsystems

### Average MLC leakage

### Primary collimation

### Jaw collimation

- Travel range
- Axis of travel
- Speed of travel

- Basic dimensional description

### Multileaf collimation

- Number of leaves
- Basic dimensional description
- Mode of operation
- Axis of travel
- Speed of travel
- Resolution
- Leaf drive mechanism
- Position indicator

0.25% (typical)

Tungsten block with rectangular, fixed aperture

1.0 cm - 5.0 cm treatment field width at isocenter  
IEC-y

Selectable fixed primary collimation widths set  
at the start of a given treatment

13.5 cm tungsten thickness

64 binary interlaced leaves (tongue and groove side profile)

10 cm leaf thickness in beam direction

Binary leaves either fully in or fully out of beam path

1 axis, longitudinal direction (IEC-y)

Binary leaf state changed within 20 msec

0.625 cm leaf widths in IEC-x direction at isocenter

Pneumatic

Continuous monitoring with interlock, checks that leaves  
open and close at correct time

## 6. Dosimetry

### SPECIFICATION

Chamber type  
Precision  
Linearity  
Variation of dose rate vs. gantry angle  
Beam-off interlocks

### PERFORMANCE

Full beam transmission, hermetically sealed dual chamber  
Within 1%  
Within 1%  
Within 2%  
Radiation System; Dose Monitoring, Patient Table Position; Gantry Angle; Computer Communications; Vault Door; Facility E-stops

## 7. CTrue™ Imaging

### SPECIFICATION


Geometry  
Dose per MVCT image (typical)  
Detector configuration  
  
Image resolution (xy)  
Slice spacing available  
Scan time  
Field of view (FOV)  
Source to detector distance  
Isocenter to detector distance  
Image noise  
Image uniformity  
Spatial resolution  
Contrast resolution  
Image reconstruction algorithm  
Image reconstruction time  
Image registration tools  
  
Application of adjustments obtained via image registration  
Frequency of correction for geometry flexion  
Image guidance mode

### PERFORMANCE


Fan-beam  
0.5 - 3 cGy depending on resolution and body thickness  
528 channels, single-row xenon ion chamber array used for image acquisition  
512 x 512 (0.78 mm pixels)  
2 mm, 4 mm, 6 mm  
Typically 2 minutes per 10 cm length at 4 mm slice spacing  
40 cm diameter  
145 cm  
60 cm  
2% - 4%  
Within 2.5%  
Nominal 0.5 lp/mm at 10% MTF  
2% density for 2 cm object (typical)  
Filtered backprojection  
Real-time; slice-by-slice at time of acquisition  
MVCT/kVCT image overlay with adjustable checkerboard; manual or automatic registration (maximization of mutual information) using bone and/or soft tissue; translations and pitch/roll/yaw determination  
Translations applied via couch, roll applied via gantry  
  
None required; rigid ring gantry platform  
Daily 3D MVCT matched with 3D kVCT




Prostate visibility

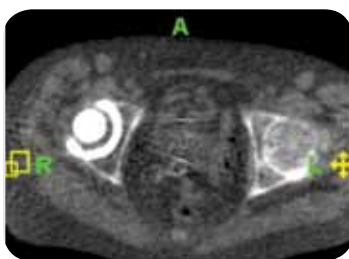


No high-Z artifact

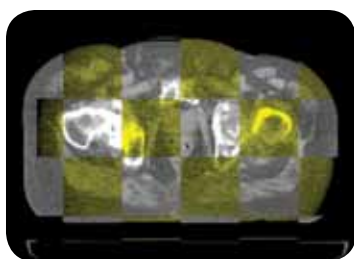

Multi-planar image overlay

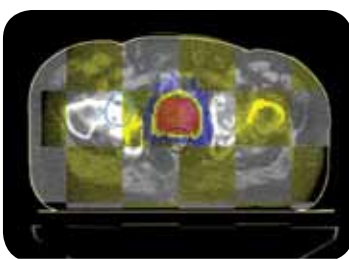

Overlay images, structures, dose

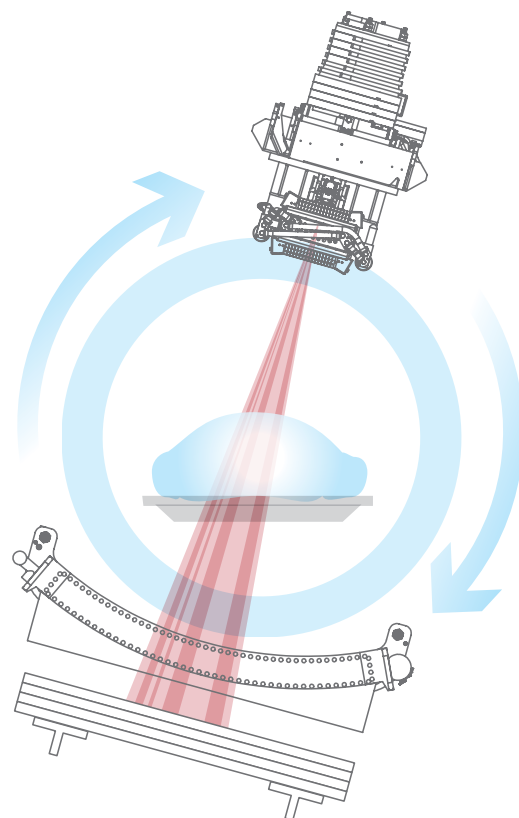

Common imaging and treatment subsystem

## 8. Safety Features

**Integrated treatment planning:** Planning and treatment systems use a common database

**Integrated quality assurance:** Patient-specific QA via a dedicated application on the planning station

**Integrated R&V:** Built-in record and verify in addition to DICOM compliant MOSAIQ™ and Aria™ interfaces

**Single database:** Common database for planning, imaging and treatment data

**No plan data transfer:** Planning data retrieved by treatment system from common database

**Common imaging and treatment beam source:** Imaging geometry consistent with treatment geometry (imaging and treatment source are the same), thus avoiding the need for repeated geometric calibration of 2 different sources

**Fully enclosed gantry:** Integrated gantry enclosures protecting the patient from beam delivery system components

**No auxiliary beam collimation:** Integrated binary MLC provides all beam collimation and modulation, removing the possibility of incorrect component installation and maintaining maximum physical patient clearance; also reduces handling of hazardous material

**No bending magnet:** No beam bending required, therefore eliminating the possibility of geometry errors resulting from redirection of the beam

**No electron mode:** The treatment system cannot be set to an unintended mode of delivery (note that helical IMRT with 6 MV x-rays can be planned to deliver treatments similar to electron mode treatments)

**Full system integration:** Level of system integration results in software for all components being designed and tested in unison, reducing concerns of disparate system interfaces and compatibility concerns across products and across companies



## STANDARD SYSTEM


### 9. TomoHelical™ Treatment Delivery Mode

The *TomoHelical* delivery mode provides IMRT and 3D CRT treatment delivery in a continuous (360°) helical pattern, using thousands of narrow beamlets, which are individually optimized to target the tumor. The *TomoHelical* mode maximizes conformity and uniformity of dose to the tumor while minimizing exposure to healthy tissue.

The user is able to create a treatment plan that defines dose goals and constraints for target and avoidance structures, the level of modulation for the plan, as well as the fractionation schedule.

During treatment delivery, the linear accelerator completes multiple 360° rotations around the patient while the couch passes through the bore of the system, initiated by a single turn of the operator console key.

Targets of up to 135 cm in length\* can be treated, with no need to reposition the patient and with no field junctioning.

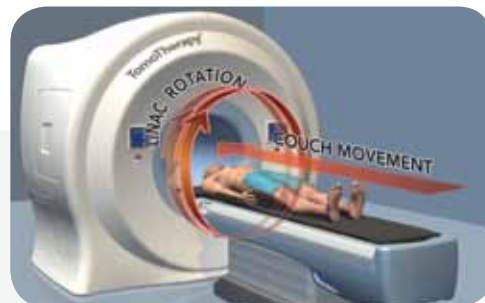

*The gantry rotates around the patient in TomoHelical delivery mode, delivering radiation in a continuous spiral pattern*


### 10. TomoDirect™ Treatment Delivery Mode

The *TomoDirect* delivery mode is a discrete angle, non-rotational delivery mode.

*TomoDirect* allows creation of treatment plans that include between 2 and 12 target-specific gantry angles. It also allows the user to define the level of modulation for the plan, including a 3D delivery mode. Treatment planning is completed rapidly due to the power of the system's computing platform.

During treatment delivery, all beams for each target are delivered sequentially with the couch passing through the bore of the system at an appropriate speed for each gantry angle. The complete treatment delivery is initiated by a single turn of the operator console key.

Targets of up to 135 cm in length\* can be treated, with no need to reposition the patient and with no field junctioning.

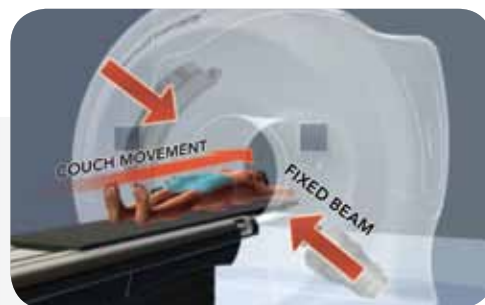

*2 – 12 discrete gantry angles can be used in TomoDirect delivery mode*

### 11. High Performance Couch

The High Performance Couch, with Medical Intelligence indexing system, provides sub-millimeter accuracy and precision in point-to-point and translational positioning.

Clinical workflow is enhanced with ergonomically designed dual Couch Control Keypads mounted to each side of the couch. The Couch Control Keypads allow motorized patient position modification in the X, Y and Z directions with simple, single-handed operation.

The custom patient couch has a high strength carbon-fiber top with an indexing system designed to accommodate immobilization systems from multiple vendors.

*See High Performance Couch brochure for additional specification information.*

### 12. Laser Positioning System

The *TomoHD* system configuration includes stationary green lasers for virtual isocenter and moveable red lasers for patient positioning and registration.

\* Typical patient set-up; Actual treatment volume lengths are variable depending upon couch height.



## 13. Computing Hardware

Includes the following items, housed in a rack enclosure (dimensions: 78"/198 cm H x 26"/66 cm W x 38"/97 cm D):

- The High Performance Optimization Cluster, which comprises a Blade Server chassis plus 14 server blades (used for treatment plan optimization and dose calculation);
- The Data Server and Storage Area Network (SAN), which stores data required to run the treatment system and deliver prescribed treatments to the patient. It also stores patient planning information (CT images, contours, dose volumes, and plans). The Data server is connected to the Optimization Server, Planning Stations, Operator Station and the Treatment Delivery System;
- The TomoGateway™ system hardware enables remote system diagnostics and monitoring by TomoTherapy Technical Support (note: additional software may be required which must be purchased separately);
- Tape drive system for database backup;
- KVM and LCD display for administration of the Optimizer/Data Server assembly;
- Uninterruptible Power Supplies (UPS) to support the complete Optimizer/Data Server assembly;
- Network switch and firewalls.

## 14. Planning Stations

Two TomoTherapy® Planning Stations allow the definition and management of treatment plans and initiation of plan optimization. Delivery Quality Assurance (DQA) tools are also integrated into the Planning Station software. The Planning Stations further include the Data Management System application software, for archiving and management of patient data. Each Planning Station is provided with a color printer capable of printing treatment plans. An Uninterruptible Power Supply, high-resolution LCD monitor, keyboard, mouse, and required cables are included with each Planning Station.

*An optional third Planning Station is available for purchase.*

## 15. Operator Station

One TomoTherapy® Operator Station, which connects to the *TomoHD* treatment system database and provides control of CTrue™ MVCT imaging, patient treatment and quality assurance delivery, and basic record & verify functionality.

The Operator Station is provided with a color printer capable of printing CTrue images and treatment data, plus an LCD monitor, keyboard, mouse, and required cables.

## 16. TomoPortal™ Remote Viewer

The TomoPortal™ Remote Viewer securely and easily provides a web-enabled link to patient information stored in the *TomoHD* treatment system. It is possible to review plan, registration and treatment data from down the hall or across a continent. Price includes server hardware and software required to run up to 2 concurrent system users.


*Options to purchase up to six additional TomoPortal Remote Viewer licenses are available.*

## 17. Standard System Interfaces

### OIS Connect™ Software

The *OIS Connect* software provides the ability to interface a *TomoHD* treatment system to a compatible Oncology Information System (OIS). The software facilitates greater integration of the *TomoHD* treatment system into the radiation oncology department, by:

- Allowing scheduling of TomoTherapy<sup>SM</sup> treatments on the OIS;
- Providing automatic capture of *TomoTherapy* procedures on the OIS;
- Aiding in charge capture and billing (where applicable);
- Aiding in integrating *TomoTherapy* treatments into patients' electronic medical records, via the OIS.

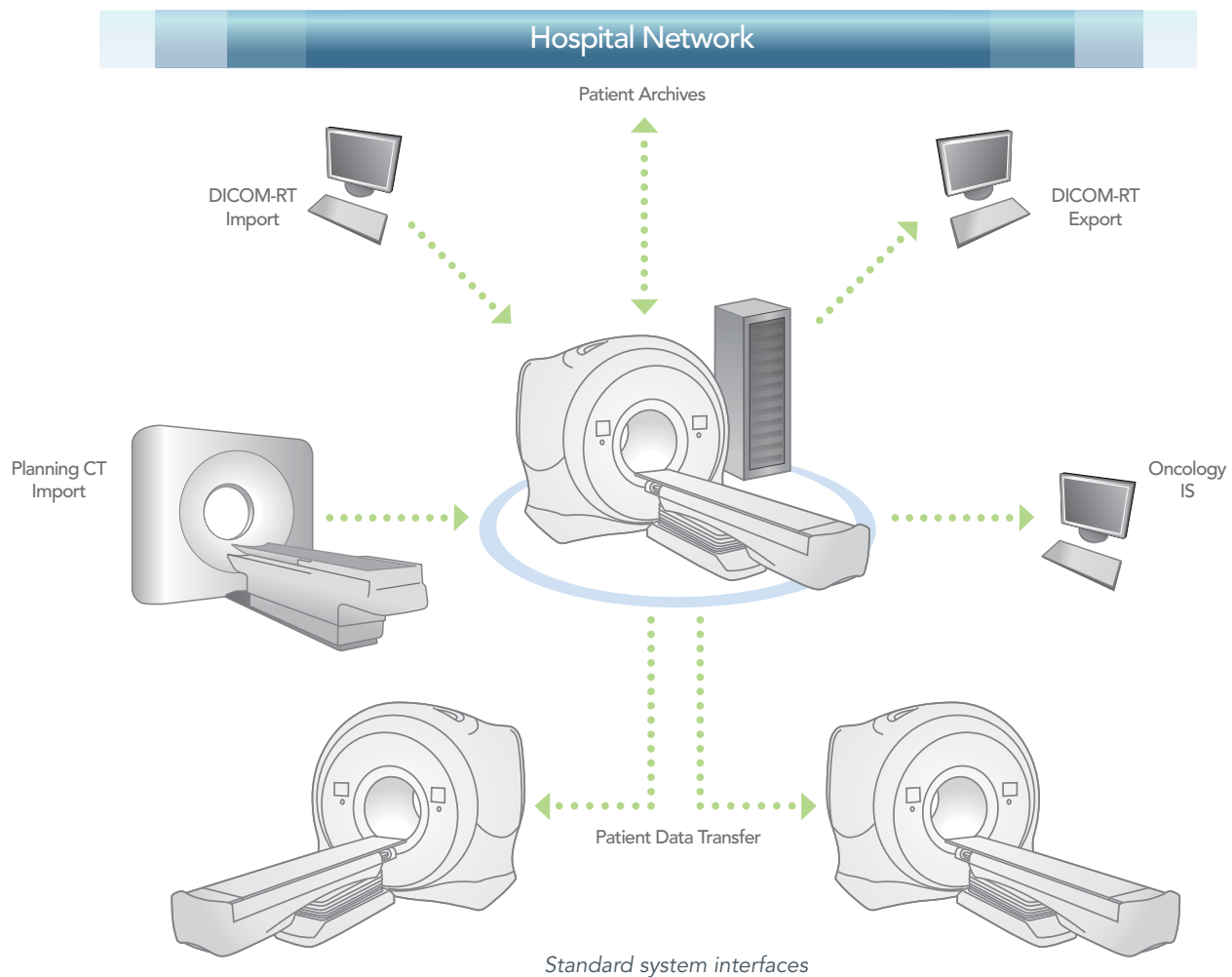

The *OIS Connect* software is based on DICOM-RT Worklist communication, as specified in DICOM Supplements 74 and 96.

### DICOM Export

The DICOM Export Data Services Package allows the following DICOM objects to be sent from the *TomoHD* treatment system to third-party systems and clinical/research databases:

- DICOM-CT Image Set
- DICOM-RT Structure Set
- DICOM-RT Dose
- DICOM-RT Plan

### Network Data Storage

This feature allows patient archives created on the *TomoHD* treatment system to be sent to, or retrieved from, a storage location outside of the treatment system network.

Network Data Storage provides configuration of the system firewall, workstation configuration and FTP system set up for data transfer via network.

### Patient Data Transfer

Patient Data Transfer allows the transfer of patient data from the *TomoHD* treatment system upon which it is installed, to another *TomoTherapy* treatment system.

Patient Data Transfer provides configuration of the system firewall, workstation configuration and FTP system set up for data transfer via network, and an external USB hard drive for physical data transfer.

## 18. Quality Assurance

### QA Package (Plan quality assurance)

The standard QA package includes quality monitoring and performance testing phantoms. Components include the following:

- One Slit Beam Virtual Water;
- One TomoTherapy® Phantom & Holder ("Cheese Phantom") with Density Plug Set;  
12 density plugs are included, representing the following materials: LN-300 Lung, LN-450 Lung, Adipose, Breast BR12, Brain, CT Solid Water, Liver, B200 Bone Mineral, CB2-30% CaCO<sub>3</sub>, CB2-50% CaCO<sub>3</sub>, Cortical Bone and Inner Bone;
- Two Calibrated Mini Ion Chambers, and approximately 59' (18 m) Interconnect Cables;
- One Calibrated CT Slice Ion Chamber, Buildup Cap and Jig.

### Tomo® Quality Assurance (TQA™) Total Package (Machine quality assurance)

The TQA application is a calendar-based productivity tool that simplifies the collection and analysis of machine performance data for the TomoHD treatment system.

The application leverages internally-generated data to provide results quickly and easily. The TQA application offers trending and reporting of many system and dosimetric parameters that allow physicists to monitor the performance of the TomoHD system. All data may be exported.

The Total Package includes the base TQA package, Enhanced Dosimetry package, Longitudinal Beam Profile package, and the LINAC Alignment package, plus the TQA Step-Wedge Phantom.

*See TQA brochure for more information.*

## 19. Installation Services

Includes:

- Pre-installation Site Planning and Project Management services;
- Installation and commissioning;
- Completion of Acceptance Test Procedure (ATP) and system handover.

## 20. Partner Marketing Portal

The Partner Marketing Portal web site is offered as an online resource to assist TomoTherapy customers in promoting their investment in cancer care. Marketing resources available via the site include:

- Printed patient and referring physician materials available for order;
- Full image library, including product and patient photography;
- Embeddable animations and video for treatment centers' websites;
- Customizable TomoTherapy-focused TV spots;
- Suggested press release template for technology introduction.



Register for access at [TomoTherapy.com/marketing](http://TomoTherapy.com/marketing)

# ***TomoTherapy®***



## UNITED STATES

### **Accuray Corporate Headquarters**

1310 Chesapeake Terrace  
Sunnyvale, CA 94089  
USA  
Tel: +1.408.716.4600  
Toll Free: 1.888.522.3740  
Fax: +1.408.716.4601  
Email: [sales@accuray.com](mailto:sales@accuray.com)

### **TomoTherapy**

1240 Deming Way  
Madison, WI 53717  
USA  
Tel: +1 608 824 2800  
Fax: +1 608 824 2996

## ASIA

### **Accuray Japan K.K. - Tokyo**

Daini Tekko Building 6F  
1-8-2 Marunouchi, Chiyoda-ku  
Tokyo 100-0005  
Japan  
Tel: +81.3.6269.9556  
Fax: +81.3.3217.0337

### **Accuray Asia Ltd. - Hong Kong**

Suites 1702 - 1704, Tower 6  
The Gateway, Harbour City  
9 Canton Road, T.S.T.  
Hong Kong  
Tel: +852.2247.8688  
Fax: +852.2175.5799

### **Chengdu Twin-Peak Accelerator Technology**

#39 Huatai Road, Longtan  
Industrial Zone,  
Section 2 East, 3rd Ring Road,  
Chengdu, Sichuan 610051  
China

## EUROPE

### **Accuray Europe**

Tour Atlantique 25e  
1 Place de la Pyramide  
92911 Paris La Défense Cedex  
France  
Tel: +33.1.5523.2021  
Fax: +33.1.5523.2039

### **TomoTherapy Belgium BVBA**

Pegasuslaan 5  
1831 Diegem  
Belgium  
Tel: +32 (0)2 400 4400  
Fax: +32 (0)2 400 4401/02